




Power, Precision, Performance

OEC 9900 Elite



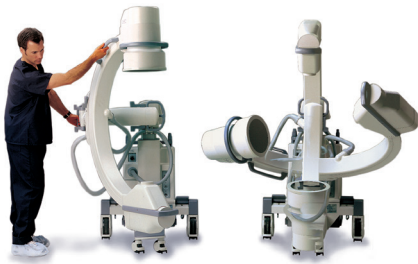

Maneuverability and Positioning

- OEC 9900 platforms have 31" of free space in between the X-ray tube and the image receptor to provide a large imaging space.
- The C-arm depth is 26" or deeper to provide a large imaging space and C-arm clearance around the patient and the imaging table.
- The C-arm models provide a minimum of 115° C-arm orbital rotation, 90° under-scan and 25° over-scan. The Super-C model increases the over-scan to a maximum of 55°.
- The system mechanics allow the user to reverse the X-ray tube and the I.I. positions and also provide C-arm under-scan and over-scan capabilities.
- The C-arm is able to rotate $\pm 180^\circ$ to allow the imaging chain to accomplish angled projections.
- The system has at least 18" of vertical C-arm travel capability to adjust the imaging chain height.
- The C-arm provides side-to-side (wig-wag) and the horizontal travel movements to allow panning during imaging.
- The system is equipped with cable pushers on the C-arm wheels so that the C-arm can more easily move over the cables.
- The system can be upgraded to provide motorized C-arm movements in both orbital and lateral directions at 9° per second.




Power to Perform

- The generator is a 60 kHz high frequency inverter.
- The output power rating of the generator is 15 kW.
- The system operates at full capacity out of a standard wall outlet on 120 Volts AC, 15 amps.

- The generator meets the following power requirements:
 - Fluoroscopic kVp range: 40 - 120 kVp
 - Fluoroscopic mA range: 0.2 - 10 mA
 - Radiographic kVp range: 50 - 120 kVp
 - Radiographic mA range: up to 75 mA
- The generator is capable of providing a boost or a high dose fluoroscopic exposure at up to 20mA.
- The generator is capable of providing pulsed fluoroscopy at 1, 2, 4, or 8 pulses per second.
- The generator is designed to provide a high mA digital spot or digital radiography mode with up to 75 mA.
- The generator is capable of pulsing at 1, 2, 4, or 8 pulses per second in normal fluoroscopy mode to reduce radiation exposure if needed.
- The generator is capable of pulsing at 30 pulses per second in Digital Cine Pulse Mode to reduce imaging lag caused by patient or C-arm movement.
- The maximum X-ray tube current in radiographic mode is 75 mA from standard 120 Volts AC, 15 amp outlet.
- The mAs range in radiographic mode is between 1 to 300 mAs.

It's All About Imaging

The Tube . . .

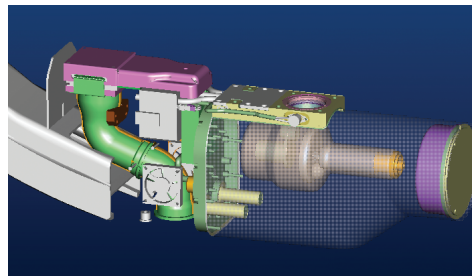

- The X-ray tube uses a rotating anode X-ray tube.
- The X-ray tube has dual focal spots sizes of 0.3 mm & 0.6 mm nominal.
- Anode heat storage capacity is 300,000 H.U.
- The anode cooling rate is 85,000 H.U./min.
- Anode temperature is computer monitored for protection. An on-screen warning is provided before the anode reaches its maximum heat storage capacity.
- The tube housing heat storage capacity is 1.6 million H.U.

- Proprietary heat pipe cooling technology and a ducted fan cooling system improve efficiency of X-ray tube heat management so the system can meet the requirements of higher power interventions.
- Protection is provided for tube overload.

The Image Intensifier . . .

- The system has a 9" tri-mode image intensifier.
- The minimum central resolution at the monitor in the 9" mode is 2.2 line pairs/mm.
- The system has a 12" tri-mode image intensifier option for vascular applications.
- The minimum central resolution at the monitor in 12" mode is 1.6 line pairs/mm.
- The image intensifier is equipped with a carbon fiber grid with 10:1 grid ratio and 150 lines/inch.

Display Details

- The system is equipped with a high resolution 1K x1K CCD camera.
- The camera gain and iris are computer controlled.
- The system is equipped with an automatic ABS tracking system to optimize image quality.
- The system is equipped with two 19" LCD high resolution monitors.
- The system has an articulating monitor arm capable of moving 22" in the horizontal direction and a 7° up /10° down degree tilt along the vertical plane. Monitors are viewable from all four sides of the workstation.
- The system provides an ambient room light sensor to automatically adjust the monitor brightness for optimum image display.
- The system has a touch screen user interface.
- The system has a reference monitor on the C-arm which displays the live fluoro image. (Not available on the motorized system.)



Image Management and More

- The system provides a simple method to input patient information into the system computer so that each image is properly annotated.
- The system allows the user to enter multiple patients into the system via a work list.
- The system provides a last image hold capability where the last image is displayed on the active monitor after the termination of an exposure.
- The system offers Dynamic Range Management (DRM)- Precision Imaging and Pre-set Imaging Profiles.
- The system is equipped with a backlit X-ray control panel that allows technologists to operate the system in a dimmed operating room.
- The system allows the user to change the image orientation on the display screen during a live exposure or using the last image hold. Those functions include image rotation, left to right and top to bottom image reversals.
- The system is equipped with an iris collimator and shutters to limit the imaging field and reduce the X-ray exposure area to improve image quality.
- The system provides an on-screen collimator position indicator which can be applied on the last image hold so that the user can adjust collimation without X-ray exposure.
- The system allows the user to change the image field of view in real-time to magnify images and improve image resolution for better visualization of small anatomical parts.
- The system provides a means to monitor accumulated exposure time during a procedure and an alarm to warn the user of long exposure time.
- The system provides means to record and display radiation exposure dose information for a procedure.
- The system automatically selects proper imaging technique, kVp and mA, during imaging, but also provides the user with the ability to override these settings manually.

- The image processing matrix is 1,000 x 1,000 x 16 bit.
- The system is capable of providing digital cine modes at either fifteen (15) or thirty (30) pulses per second.
- The system provides noise reduction control to improve image quality.
- Real-time automatic brightness and contrast control are provided to optimize the displayed image.
- The system provides a real-time and post processing edge enhancement capability to improve bones and vascular details.
- An electronic zoom function is provided to magnify a captured image for improved visualization.
- The system provides an automatic save function so that the last image taken is automatically saved to the hard disk.
- The system provides means to put markers or comments on an image.
- The system provides a simple method to retrieve stored images (static and dynamic), display a collage of stored images, and allow the user to review stored images at a later time.
- The system is capable of saving 63 - 1000 static images to the internal hard disk, depending on the product configuration. Stored images can be retrieved at a later time.
- The system has an integrated CD/DVD burner and allows saving in DICOM or OEC format.
- The system has a USB 2.0 mass storage write-only capability to save images in .jpg/.bmp/.ave file formats. USB storage devices must be unencrypted and unprotected.
- The system is capable of storing angiographic runs (dynamic images) to the internal hard disk at various frame rates. The maximum frame rate is 30 frames per second.
- The system provides dynamic image storage of up to 10 minutes at 30 frames per second.
- The system provides an automatic playback capability so images can be played back immediately after the completion of a dynamic imaging run.
- The system allows the user to review a dynamic imaging run frame-by-frame and to set cues in the run.
- The system provides digital subtraction and roadmap imaging modes for vascular applications. These imaging modes can be accessed and operated by the user.
- The system provides pixel shift, landmark, peak opacification and re-masking capabilities in the vascular mode.
- The system provides integrated high quality hardcopy features to a thermal printer, laser quality film printer, or interface to a laser camera.
- The system provides a DICOM 3.0 interface capability that can be connected to the hospital's network system to transfer images for image archive and print purposes.
- The system's DICOM 3.0 interface provides storage class, print class and query/retrieve work list features.
- The system allows the option for wireless DICOM which utilizes a Wireless LAN 802.11 b/g/n, 10/100 MB device, runs a custom embedded Windows implementation, and currently supports WEP-64, WEP-128, WPA-PSK (TKIP), WPA2-PSK (AES-CCMP) and PEAP with EAP-MS-CHAP V2 security methods.

Training Advantages

- The purchase of this product includes in-house training of staff technologists by an applications specialist (certified technologist) from the manufacturer.



- The training includes continuing education credits for staff technologists who complete the course work.



©2018 General Electric Company – All rights reserved.

General Electric Company reserves the right to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation. Contact your GE Representative for the most current information.

GE, GE Monogram and OEC 9900 Elite are trademarks of General Electric Company. GE OEC Medical Systems, Inc., doing business as GE Healthcare.